

Hanse-Bathing-System



Relaxed Bathing · Safe Transfer · Easy Handling

EXPERIENCE THE FUNCTION – DISCOVER THE SPECIAL

hanse
MEDIZINTECHNIK

Hanse-Bathing-System



A new bathing experience

In modern patient care, bathing is not just a matter of cleanliness for those dependent on care but also provides relaxation and therapy. Hanse Bathing Systems were designed with these needs in mind.

Unlike conventional tubs, the Hanse step tubs allow patients to bathe in a safe, relaxed and comfortable seated position, rather than while lying back. On top of this, an adjustable whirlpool function provides a stimulating water current. Altogether this makes for a gentle and soothing bathing experience.

Excellent hygiene


The Hanse Bathing System's tub features a microscopically-smooth synthetic surface, without any joints or corners that are susceptible to dirt. The use of stainless steel provides a good appearance and excellent hygiene. Thanks to the fully-automatic disinfection system, the whole system can be automatically cleaned and disinfected at the push of a button. The precise dosage of disinfectant makes it economical in use and maintaining the highest standards of hygiene.

Hanse gets the bathing experience moving

The unmistakable whirlpool effect in Hanse baths is provided by their twin-stream system, allowing the warm water to flow all around the body offering an effortless bathing experience while at the same time being both invigorating and gentle on the circulation. The whirlpool function is adjustable and can be further enhanced with air bubbles, providing an even more intensive stimulation of the patient's circulation.

The combination of the special seat and water channelling system ensures that the water flows all around the body, cleansing thoroughly even in sensitive intimate areas. Adding high-quality bath preparations enhances the whirlpool's cleansing, revitalising and skin-care effect even further.

The whirlpool system also has a connection for an underwater massage hose for applications with suitably trained care personnel.

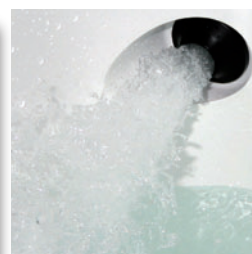

Hanse-Bathing-System

Simple handling

The entire transport and bathing process is back-friendly for the carer and can be operated by just one member of the nursing staff.

The integrated lift can be raised and lowered smoothly, even with heavy patients, and its simple operation helps make life easier for nursing staff.

Thanks to the operation by foot pedal, the nursing staff can focus on the patient.



focus on the patient



transfer to the bathroom



clearly arranged fittings



safe hoisting and lowering



ergonomical drying

Safety at all times

Thanks to the anatomically-formed seat, padded pelvic belt and padded armrests that can be raised individually, the patient can feel secure even during transport to and from the bath.

Because the height of the seat is matched to that of the hospital bed, this ensures that the patient can be transferred simply and safely. The seat is also ideal for showering and positioning over a toilet.

Models

You can choose from 4 models

Whether in retirement homes, disabled care or as an wellness bathing tub in a snoezelen room – the Hanse Bathing System complies highest requirements.

For a nice spot of colour in the bathroom the tub-mantle is available in different colours and also with effect coating.

Space-saving design

A further advantage of the bathing system is its space-saving design, allowing efficient use even in spaces as small as 1,70m x 2,50 m upwards, something that would be unimaginable with a conventional bathing system.



Type 125

- Step tub
- Integrated patient hoist with pivoted body-shaped seat
- Longer shape for mobilization exercises



Type 95

- Step tub
- Integrated patient hoist with pivoted body-shaped seat
- Even for smallest bathroom from 1,70m X 2,50m



Type 125 T

- Step tub with body-shaped seat
- Electromotive operating side door
- Longer shape for mobilization exercises



Type 125 TL

- Step tub
- Electromotive operating side door
- Integrated patient hoist with pivoted body-shaped seat
- Longer shape for mobilization exercises

Technical Data:

Max width:	720 mm
Max length Type 95:	1.480 mm
Type 125 (T/TL):	1.780 mm
High (top of the tub):	930-960 mm
Max Load:	150 kg

Detailed information see system description.

Accessories:

- twin whirl action unit
- automatic disinfection device
- host for underwater massage
- high-quality shower attachment set
- transport chair · bucket for transport chair · head rest
- seat cushion with hygiene opening
- tub mantle in RAL-colour / effect-coating